Desiccators

Model : SDC-30(U),45(U)

Manual No. : D5011L002 Version :4.1

1. Characteristics

(1) Multi-stackable Design

The SDC models come in two capacities to suit the user's needs with the internal volume ranging from 35 liter and 50 liter.

With the same foot print, different models can be stacked together in any convenient fashion allowing for lab space saving.

(2) Efficient Use of the Internal Space

Compared to round desiccators, cubic desiccators provide larger usable volumes with the same footprint. In addition, the internal space can be more efficiently utilized by using the standard perforated shelves designed for efficient airflow.

Two perforated shelf or three perforated drawers are provided and, if needed, up to six shelves(option) can be used for efficient storage of the samples. (Refer to Product Checklist or Specifications for details.)

The level of shelves or drawers can be adjusted by built-in rails.

(3) Easy Humidity Level Checking

With an built-in Thermo-hygrometer installed on the front door, the humidity level checking can be done conveniently.

(4) Integral Molding of the Polycarbonate Body

Unlike conventional vacuum desiccators, the polycarbonate cubic body (except the door part) of this unit is integrally molded as a single piece for maximum sealing performance.

Shock-resistant, high-tensile strength, and transparent polycarbonate body provides extreme durability as well as 360° visual observation from outside.

(5) Desiccant Tray

This unit comes standard with a desiccant tray.

(6) UV block effect (only for U models)


SDC-30U, 45U models minimize damages or discoloration of light-sensitive samples.

2. Precautions

CAUTION

- Note that this unit cannot be autoclaved.
- When putting high-temperature samples inside the unit, cool them thoroughly in advance or make sure that they are not in direct contact with the interior of the unit to avoid interior damage.
- Do not move this unit when it contains sample inside.
- Do not use chlorine bleach, ammonia-based cleaners, abrasives, organic solvents, or metal scouring pads when cleaning. Use a soft cloth all the time.

3. Layout

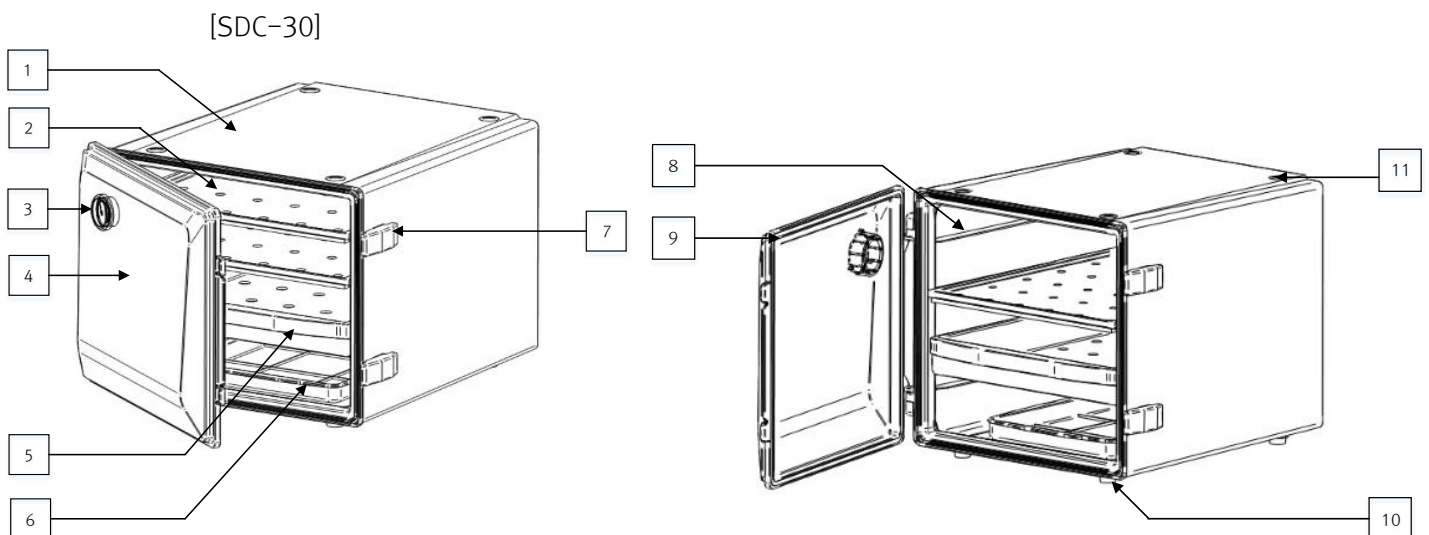- (1) Body
- (2) Perforated shelf (option)
- (3) Thermo-hygrometer
- (4) Door
- (5) Perforated drawer
- (6) Desiccant tray

- (7) Door latch
- (8) Rail
- (9) Door seal
- (10) Non-slip feet
- (11) Stacking groove

4. Thermo-hygrometer Separation and Installation


- (1) When replacing the battery, separate the thermo-hygrometer attached to the case from the product by pulling it forward. Then remove the thermo-hygrometer from the case.

- (2) Open the battery cover on the back and replace the battery.
The battery type is LR44 that can be used continuously for about 4 months.

- (3) After replacing the battery, install the case and thermo-hygrometer again in reverse order. Then push thermo-hygrometer into the place at the back side of the door following arrow direction as below figure.

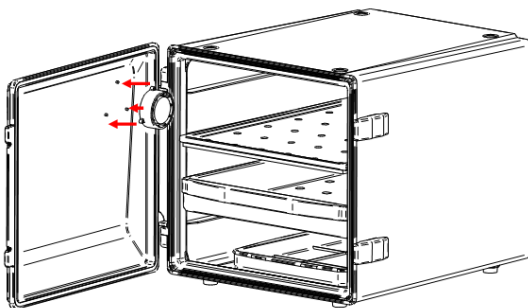- (4) Close the door and check thermo-hygrometer is installed in the right place.


NOTICE

- The thermo-hygrometer on the front can be used for reference only. If you need accurate humidity data, use an appropriate instrument or sensor.
-

5. Shelves and Drawers Placement and Separation

SDC models come with perforated drawers and a desiccant tray. Follow below steps to adjust perforated shelf or drawer. If needed, up to four (SDC-30) to six shelves (SDC-45) can be used for efficient storage of the samples. (Refer to Product Checklist or Specifications for details.)

- (1) Pull a shelf or a drawer to take out.
- (2) When the bump (refer to below circle) of the shelf or drawer is stuck in the rail, pull it out hard to take out.

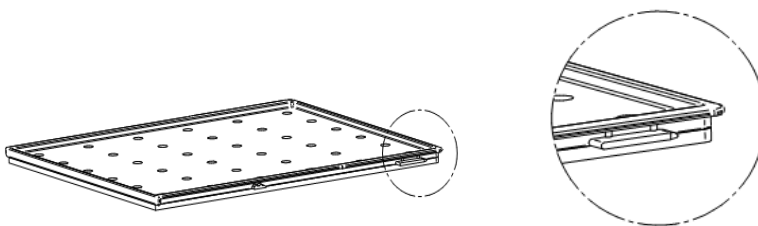- (3) Place a shelf or a drawer into any rail by following above steps (1)~(2).

NOTICE

- For user safety, perforated shelves and drawers are designed not to be taken out at once.
-

6. Specification

Model		SDC-30	SDC-45
Internal volume (L/cu ft)		35/1.2	50/1.8
Thermo-hygrometer		Digital	
Material	Body, Perforated shelves, drawer	Polycarbonate	
	desiccant tray	Polypropylene	
	Door latch	Stainless steel	
	Foot	Silicone rubber	
Shelves	Number of standard / Max.	Optional/4	Optional/6
	Max. load per shelf (kg/lbs)	5/11.0	
Drawers	Number of standard / Max.	2/3	3/5
	Max. load per drawer (kg/lbs)	5/11.0	
Dimensions	Overall (W x D x H, mm/inch)	360 x 413 x 326/ 14.2 x 16.3 x 12.8	360 x 413 x 446/ 14.2 x 16.3 x 17.6
	Internal (W x D x H, mm/inch)	322 x 367 x 302/ 12.7 x 14.4 x 11.9	322 x 367 x 422/ 12.7 x 14.4 x 16.6
	Perforated shelves (option) (W x D x H, mm/inch)	279 x 354 x 3/ 11.0 x 13.9 x 0.1	
	Drawer (W x D x H, mm)	279 x 346 x 45/ 11.0 x 13.6 x 1.8	
	Desiccant tray (W x D x H, mm/inch)	218 x 268 x 26/ 8.6 x 10.6 x 1.0	
Net weight (kg/lbs)		4.2/9.3	5.0/11.0

Model	SDC-30U	SDC-45U
Color	Amber	
UV Protection (Range)	UV A 99.8% (315~400nm)	

※ Above specifications can be changed without any prior notice.

7. Unpacking and Installation

(1) Checking Shipment Damage and Unpacking

Upon receiving the instrument, check to ensure that no damage has occurred during shipment. It is important that any damage that occurred during shipment is detected before unpacking. If you find such damage, the carrier must be notified immediately.

(2) Checking the Contents of the Package

After unpacking, check to ensure that all the parts and accessories described below are included in the package. If not, contact your dealer immediately.

Product Check list

Item	Quantity
Main Body	1
Thermo-hygrometer	1
Door seal	1
Desiccant tray	1
Perforated drawer	SDC-30(U) : 2 SDC-45(U) : 3
User manual	1

8. Warranty

1) Warranty period

- (1) If the product fails during normal and proper use, the warranty period for manufacturing liability shall be two years from the date of delivery.
- (2) If you request a repair, please check the items below and let us know so that you can get a quicker and more accurate repair.

- | |
|---|
| <ul style="list-style-type: none">• Date of purchase:• Customer name / address / phone number / E-mail:• Fault status• Model• Serial number |
|---|

2) Technical Service Contact Points

- Technical Services
- Address : 153 (Youngsan-dong), Techno 2-ro Yuseong-gu, Daejeon-city, 34025, Republic of Korea
- Tel : +82 (0)2 -2627-3824
- Website : <https://www.JeioTech.com/eng/>
- E-mail : overseas@JeioTech.com

3) Certificate of Product Warranty

This is a guide to the product quality warranty and service terms provided by JEIOTECH Co., Ltd.

All products are guaranteed to operate normally when installed under proper conditions, as specified in the user manual for each product, and used according to the intended purposes and permissible operating conditions.

The warranty period for product quality starts from the delivery date and is indicated in the catalog provided for each product or in the purchase order agreed upon with the customer. If a product malfunctions within the warranty period, JEIOTECH Co., Ltd. provides free repair services. However, among the transportation costs incurred during the repair process, any special costs resulting from the customer's site conditions (e.g., door or wall modifications, use of special handling equipment) are not covered under the warranty. Furthermore, JEIOTECH's liability under this warranty is limited to replacing the product or refunding the original purchase price at JEIOTECH's sole discretion if reasonable attempts to repair or replace the product fail.

Parts replaced during free repairs within the warranty period become the property of JEIOTECH Co., Ltd. Additionally, the warranty period is not extended due to free repairs. JEIOTECH is not responsible for any losses incurred by the customer due to delays in repairs caused by issues such as part shortages or shipping delays.

The following cases will be handled as paid repairs, even within the warranty period:

1. Non-defective cases

- * Requests to resolve issues arising from improper operating environments or abnormal installations.
- * Reinstallation requests or connections with external equipment other than the initial installation.
- * Requests for usage instructions or simple adjustments that do not require disassembly.
- * Inspections of non-defective products.
- * Cleaning or debris removal service requests.

2. Failures caused by user errors:

- * Malfunctions due to usage not in accordance with the user manual.
- * Issues resulting from improper operating environments or abnormal installations.
- * Malfunctions caused by incompatible chemical usage.
- * Malfunctions caused by irregular or improper electrical or utility supplies.
- * Damage due to product movement, storage, external impact, or dropping by the customer.
- * Malfunctions resulting from unauthorized modifications or repairs without prior written approval from JEIOTECH Co., Ltd.
- * Malfunctions caused by using consumables or optional components not specified by JEIOTECH Co., Ltd.
- * Repairs performed by personnel other than JEIOTECH or its authorized service agents.

3. Other cases

- * When the product model or serial number has been damaged, removed, manipulated, or altered.
 - * Malfunctions caused by fire, flooding, or other natural disasters.
 - * Failure due to consumable parts reaching the end of their lifespan (e.g., lamps, filters, fuses, gaskets, packings, hoses, heat transfer fluids).
 - * Natural wear and tear (e.g., discoloration, screen printing degradation).
 - * Damage or malfunction of product accessories (e.g., power cables, data cables, glass components) or options.
 - * Malfunctions caused by errors in computers, computer components, or connected network systems.
- Before JEIOTECH performs free or paid repairs or accepts the return of a product, the customer must ensure the product has been decontaminated of chemical and biological hazards.
 - Replacement parts used during repairs may differ from the original but will function correctly and be equivalent in performance.
 - JEIOTECH does not bear any responsibility for damage to stored samples, data loss, or any other direct or indirect damages caused by product malfunctions, regardless of the warranty period.
 - JEIOTECH does not guarantee issues related to computers, computer components, or network systems connected to the product, regardless of the warranty period.
 - This warranty is the sole and exclusive warranty provided. Any implied warranties imposed by applicable laws are limited to the warranty period stated herein. This warranty cannot be modified by any verbal advice, statements, or documents on behalf of JEIOTECH unless explicitly labeled as a "Warranty Amendment" and signed by an authorized representative of JEIOTECH.
 - This warranty applies only to products purchased for direct use by the customer. It does not cover products purchased for resale purposes without prior written approval from JEIOTECH.
 - This warranty does not apply to consumables and is non-transferable.

- ✧ The contents of this manual can be changed or upgraded without prior notice.
- ✧ The copyright of this manual is reserved by Jeio Tech.

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6. Specification

Model		SDC-30	SDC-45
Internal volume (L)		35	50
Thermo-hygrometer		Digital	
Material	Body, Perforated shelves, drawer	Polycarbonate	
	desiccant tray	Polypropylene	
	Door latch	Stainless steel	
	Foot	Silicone rubber	
Shelves	Number of standard / Max.	선택 사양/4	선택 사양/6
	Max. load per shelf (kg)	5	
Drawers	Number of standard / Max.	2/3	3/5
	Max. load per drawer (kg)	5	
Dimensions	Overall (W x D x H, mm)	360 x 413 x 326	360 x 413 x 446
	Internal (W x D x H, mm)	322 x 367 x 302	322 x 367 x 422
	Perforated shelves (선택) (W x D x H, mm)	279 x 354 x 3	
	Drawer (W x D x H, mm)	279 x 346 x 45	
	Desiccant tray (W x D x H, mm)	218 x 268 x 26	
Net weight (kg)		4.2	5.0

Model	SDC-30U	SDC-45U
Color	Amber	
UV Protection (Range)	UV A 99.8% 이상 (315~400nm)	

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