# Safety Waste Disposal System

Models: SF-14C, SF-14P, WT-20C, WT-20P

Manual No. : M6001L002 Version :3.0

## 1. Characteristics

(1) Minimization of Steam Emissions
The hinge-type lid minimizes the vapor emission of the waste fluids in the laboratory. It protects the health of laboratory workers by preventing inhalation of harmful gases, and reduces the risk of fire due to the flammable gas leaks.

(2) Dedicated Carbon Filter (Option)
Dedicated carbon filter is a system that removes odors and harmful gases to protect the health of the experimenter from harmful substances. (This filter can be easily replaced without removing any of the existing attachments.)


(3) Excellent Chemical Resistance
To enhance the chemical resistance, general-purpose safety funnel (SF-14P), electrostatic conductive model (SF-14C) are made of polypropylene (PP). The white-colored general waste container (WT-20P) model and black-colored electrostatic conductive chemical waste container (WT-20C) are made of high-density polyethylene (HDPE) for outstanding chemical resistance.

(4) Notifying the Waste Liquid Level (SF model)
When the safety funnel is filled with 70% (approximately 14 liters) or 80% (approximately 16 liters) of 20 liters of capacity, the hazard level indication function is activated for the users to manage it safely. Initially, the disposal system is equipped with an 80-percent display settings, but the users can change it to the 70-percent displaying settings if needed.

(5) Safe electrostatic conductive Features (SF-14C, WT-20C models)
The black-colored model is made of an electrostatic conductive materials with high-performance plastics, which prevent electrostatic explosions. The disposal system is also equipped with an earth cable.

(6) Easy cleaning
For better cleaning convenience, we separated the safety funnel, waste container, and the sieve inside the safety funnel.

## 2. Structure
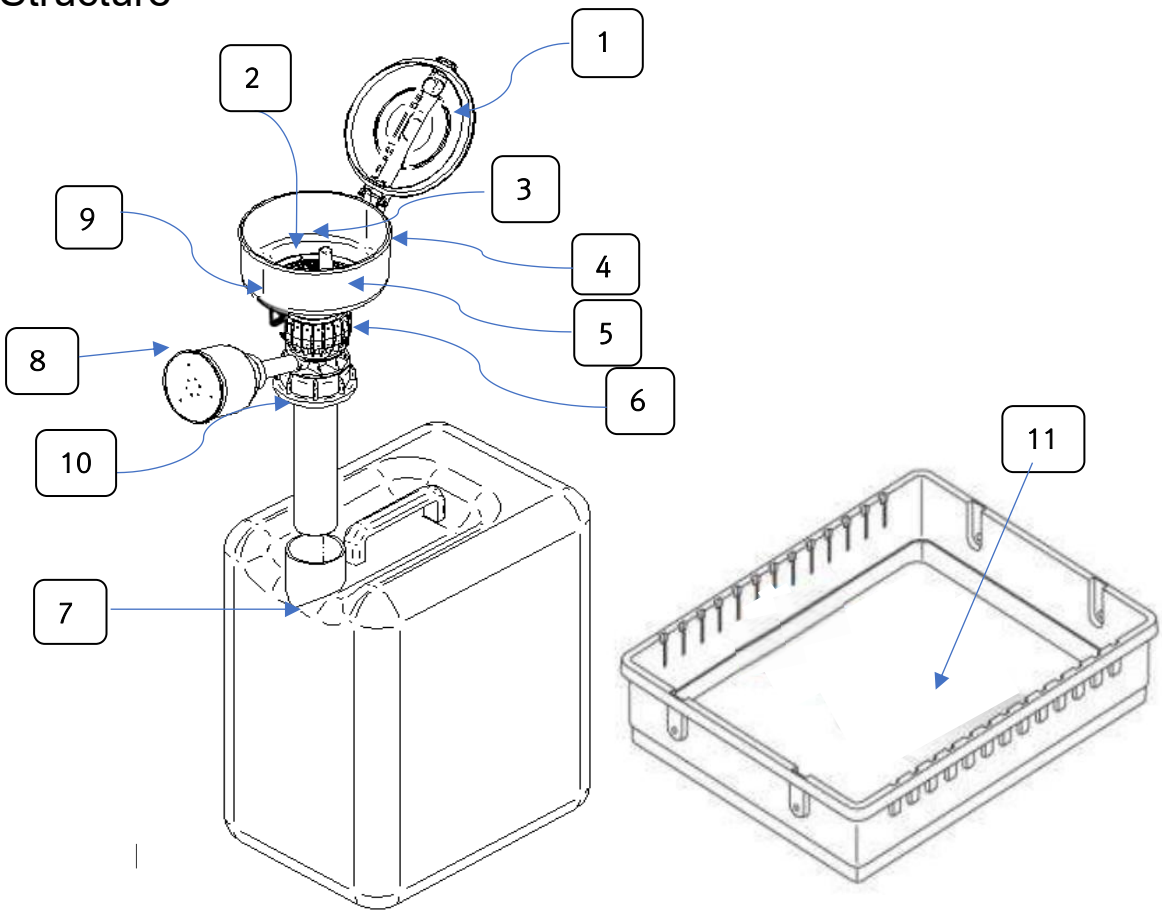
| Index Number | Name | Remark |
|---|---|---|
| **2.1** | Safety funnel cover | Hinge-type system |
| **2.2** | Sieve | Removable and washable |
| **2.3** | Level indicator | 70% or 80% level indication of the waste container |
| **2.4** | Safety funnel | Selectable between general type and electrostatic conductive model |
| **2.5** | Adapter | |
| **2.6** | Chemical waste container cap | Round perforated cap shape |
| **2.7** | Chemical waste container | Selectable between general type and electrostatic conductive model |
| **2.8** | Carbon filter | *Protection from odors and harmful substances (Option)<br>*More than 50g weight of activated carbon except for case weight |
| **2.9** | Earth cable connection | Only SF-14C model is equipped with earth cable |
| **2.10** | Level bar guide | |
| **2.11** | Tray | For storage of spillage (Option) |

©2025. All Rights Reserved. Jeio Tech Co., Ltd.

# 3. How to Use and Cautions

## 3.1. Installation and removal of the safety funnel and carbon filter

(1) When installing a safety funnel, push the container's cap (2.6) to the waste container and turn it clockwise. When removing it, turn the cap counterclockwise. (Refer to 2. Structure)

(2) When installing the carbon filter (2.8), turn one of the adapter caps installed on both sides of the adapter (2.5) counterclockwise and remove it, then slide the carbon filter (2.8) to the hole where the cap was located, and install it clockwise to the fullest extent. When removing it, turn it counterclockwise. (Refer to 2. Structure)


## 3.2. How to use safety funnel components

### (1) Sieve

The sieve (2.2) is removable. It's designed to filter magnetic bars or paper towels. If necessary, use it separately from the main body. (Refer to 2. Structure)

### (2) Replacing the Level Indicator and the chemical waste container

The Level Indicator (2.3) displays the hazard level. If the red buoy in the center is raised, please replace the waste container. There is a risk notifying sticker that reminds you to replace the waste container immediately in case the red level indicator comes up on the top of the safety funnel.

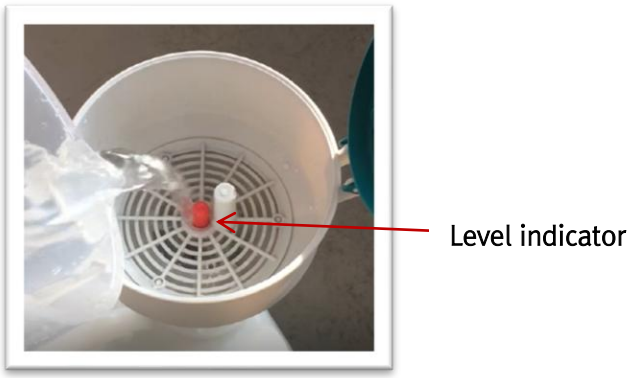**[The alarming function of the waste liquid level]**


**[Risk Notifying Sticker]**

(3) 70% level indicator

The built-in level indicator (2.3) rises when 80% of the 20 liters waste container is filled. We provide another type of the indicator that operates when 70% of the waste container is filled. Therefore replace it if necessary as shown below.

① Remove the sieve first.
② Turn the level indicator (red part) counterclockwise to pull it out.
③ Remove the level bar and separate the floater from the safety funnel where the level indicator was attached.
④ Turn the floater counterclockwise.
⑤ Replace it with a 70% level bar of 230 mm length and assemble it in a reverse order.

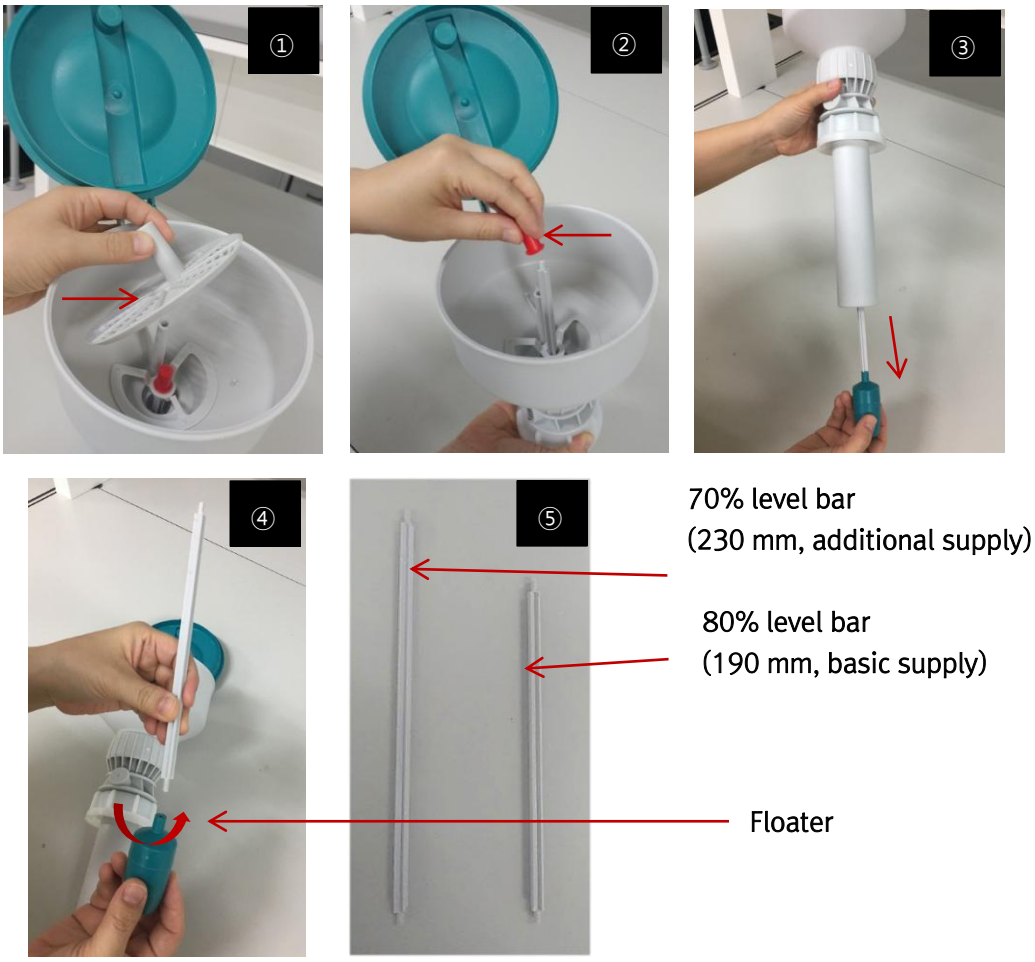※ The 70% and 80% level bar are measured based on the JEIOTECH waste container dimensions (228W x 300D x 402H, mm), and the amount of filling can be different when using a third-party waste container with different heights.


©2025. All Rights Reserved. Jeio Tech Co., Ltd.

(4) Grounding an electrostatic conductive model

When using the black electrostatic conductive model (SF-14C, WT-20C), connect the ground to prevent static fires.


- When connecting earth cable to drawers, etc.

  Use the earth cable that is built into the Earth Cable Connection (2.9).
  (Refer to 2. Structure)

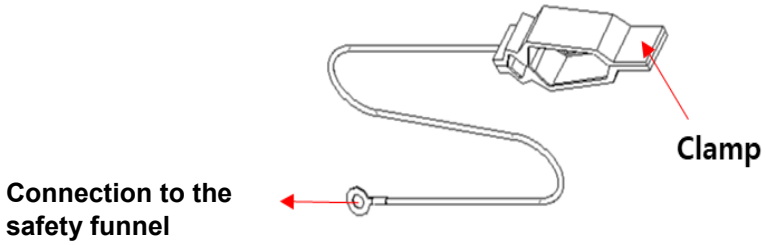**[Earth Cable]**

- Procedure for connecting the earth cable as following.
  ① Insert the enclosed screw with the earth cable into the ground connection of the safety funnel.
  ② Insert the ring-shaped end of the earth cable to the screw.
  ③ Screw the enclosed nut clockwise tightly.
  ④ Plug the clamp, which is the opposite end of the earth cable, into a suitable grounding position (such as electrostatic conductive drawers).

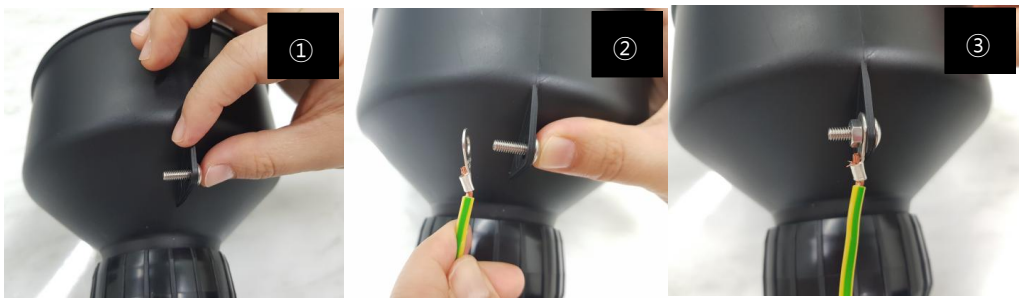- When you want to connect the earth cable built into the safety funnel directly to the outlet.

  In this case, the earth plug set must be selected as an option.
  The earth plug set provides an earth plug and an earth plug cable.


**[the front of the earth plug]**

**[the rear of the earth plug]**


**Ring-shaped end**
**: connection**
**to the clamp**

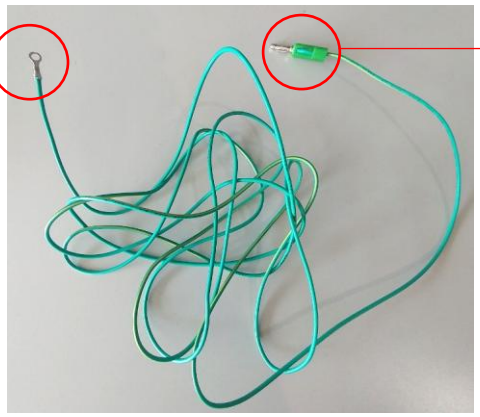**Banana plug end**
**: connection**
**to the earth plug**

[The earth plug cable]

- Procedure for connecting the earth plug cable is as follows.

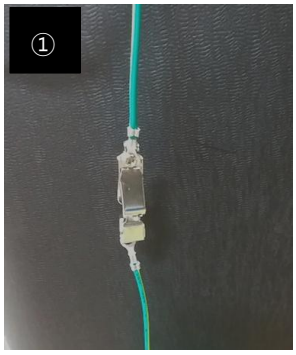① Connect to the ring-shaped end to the clamp.
(You can also connect the ring-shaped end of the earth plug cable to the safety funnel directly after removing the earth cable that already equipped to the safety funnel)

② Plug the earth plug into an outlet or an extension cord.

③ Plug the banana plug end into the earth plug.

[The earth plug set connected to the safety waste disposal system]


©2025. All Rights Reserved. Jeio Tech Co., Ltd.

(5) Tray

The optional tray (2.11) is designed to store and process waste safely in the event of a waste fluid spill. Put the waste container on top of the tray (2.11) and use it. It is highly recommended that only one waste container be stored in one tray. If two waste containers should be stored in one tray, only the same type of the waste containers must be stored together. The mixture of spillages may cause a reaction, resulting in product damage or fire.

(6) Replacing carbon filters

It is recommended to replace the carbon filter every six months. For your convenience, we provided a filter change timer to let you know when to replace the carbon filter. As shown the figures below, press the button when opening the carbon filter, then the red timer tells you when to replace the filter. You can attach it to the top, side, or anywhere you want.

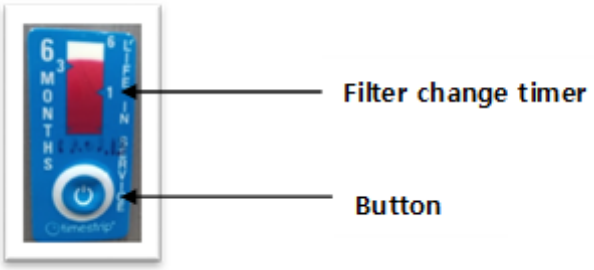**[Filter Change Timer]**

**⚠ CAUTION**

- The level indicator indicates the level when the 80% of the waste fluid is filled, so be sure to replace the waste container when the red buoy is raised.
- Connect an earth cable to the earth cable connection of the electrostatic conductive model to prevent electrostatic explosion in advance.
- Do not store two separate waste containers containing non-mixable waste fluid in one tray. Mixing may cause a reaction and cause an accident.
- It is recommended to replace the carbon filters every six months.

(7) Perforating the cap of a waste container

The waste container's cap should be perforated in a circular shape and then be installed to the safety funnel.


• If you purchase a safety funnel and a waste container together, the company will perforate the cap of the waste container in a circular shape and attach it to the safety funnel before we send it to you.

• When you buy safety funnel only, the company will send you the round-perforated waste container cap. In case, however, it does not comply with the existing waste container, please send the existing waste container's cap to the JEIOTECH dealer, then we will perforate the cap that is a circular way for you.

• If you would like to perforate the cap of the waste container by yourself, please use a hole cutter to perforate it into Ø34 mm. (Please see the figures below)


[The perforated waste container cap]

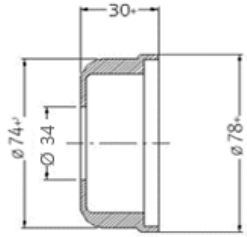[The drawings of perforated waste container]

- When you install the perforated cap,
  ① Remove the level bar guide (2.10) by turning the level bar guide counterclockwise.
  ② Remove the cap of the waste container with the nut cap by turning it counterclockwise.
  ③ Reinstall the perforated cap and then re-mounting the level bar guide (2.10) and the nut cap by turning them counter clockwise. (See the picture below)

※ Please make sure that the nut cap between the safety funnel and the cap of the waste container does not move.

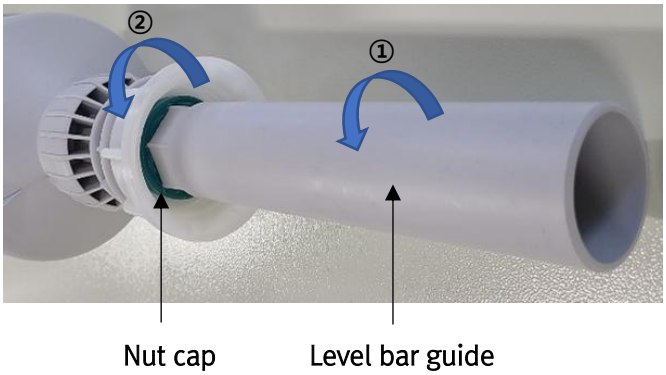[Nut cap and the level bar guide]

©2025. All Rights Reserved. Jeio Tech Co., Ltd.


## 3.3. Storage of the chemical waste container equipped with safety funnel

(1) This product is designed for safe storage in flammable cabinets of type 90 that are manufactured by JEIOTECH. It is storable in flammable cabinets with safety funnel equipped (SC3-90F-0812D2D-C, SC3-90F-0806D-C model) (see the picture below).

**[The waste container stored in the flammable cabinet]**

(2) Select the proper location by considering the size of the products. (See figure below)

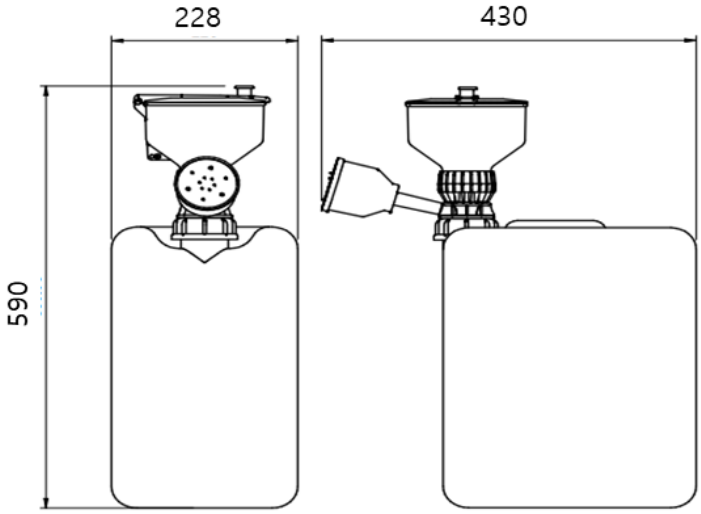**[The size of the waste container installed in the safety funnel]**

**⚠WARNING**

• When collecting the flammable substances in the waste containers, make sure to keep them in a flammable cabinet.


• To prevent static electricity, connect the earth cable before using it and use a safety funnel and a waste container of electrostatic conductive models.

## 3.4. Separation and collection waste liquid

(1) Collect the chemicals separately after checking the characteristics of them that should not be mixed.

(2) Waste containers collected from each laboratory may be sent to a larger container before they are sent to the final disposal site. The incorrect identification labels can lead to toxic vapor generation, formation of hazardous compounds, and further fire or explosion. Therefore, it is important for the users to be very careful to handle the mixture of substances that cannot coexist.

(3) This product is designed to attach the safety labels on the funnel lid so that only the chemicals you specify can be collected separately. Attach the label for identification enclosed with the product to the specified location of the safety funnel lid (see figure below).

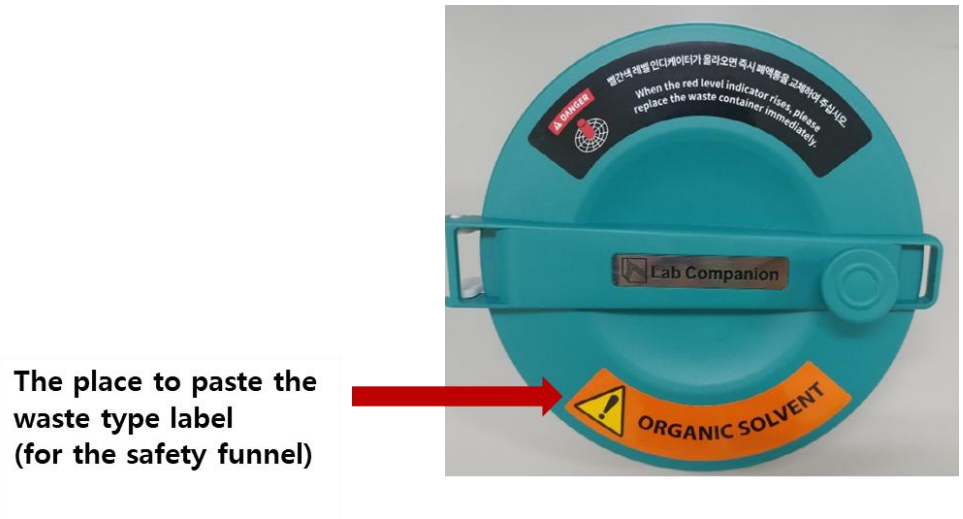[The place to paste the chemical waste identification label]

(4) There are four different colors (Persian green, Yellow, Red, Blue) of the safety funnel lid for easy

©2025. All Rights Reserved. Jeio Tech Co., Ltd.

identification of the chemicals collected. Please select one of them on your order. (Refer to 4. Special)

| Chemical waste identification labels for safety funnel | Remark |
|---|---|
|  | Collect only organic solvent waste type. |
|  | Collect only flammable waste type. |
|  | Collect only acid waste type. |
|  | Collect only alkali waste type. |
|  | Collect only halogenous waste type. |
|  | Collect only corrosive material waste type. |
|  | Collect only waste type that you have specified. |

[The type and description of the labels enclosed with the product]

(5) Labels attached to safety funnel can also be attached to waste containers to separate and collect only specified waste liquid. There are two sheets of the waste container label. Please attach the chemical waste identification label on the specified location as shown below.

**The place to paste the waste type labels (for the waste liquid bottles)**

**[The place to paste the chemical waste identification label]**

[Chemical waste container label]

폐액 식별 라벨
Chemical waste identification label

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
**HAZARDOUS WASTE HANDLE WITH CARE !!**

1. When transferring the chemical waste containers to a designated area, their caps should be tightly covered to prevent leakage.
2. The chemical waste container cannot be collected if they are damaged, the cap is not covered, or not attached with the chemical waste identification label.
3. No solid waste other than the waste liquid is allowed to be collected in the chemical waste containers.
4. Transfer the waste liquid in an iron container to the chemical waste container. Dispose of the iron container empty.
5. Always be prepared for safety accidents that may occur during the storage and transportation of experimental waste liquid.

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©2025. All Rights Reserved. Jeio Tech Co., Ltd.

[Chemical waste identification label (for chemical waste container)]
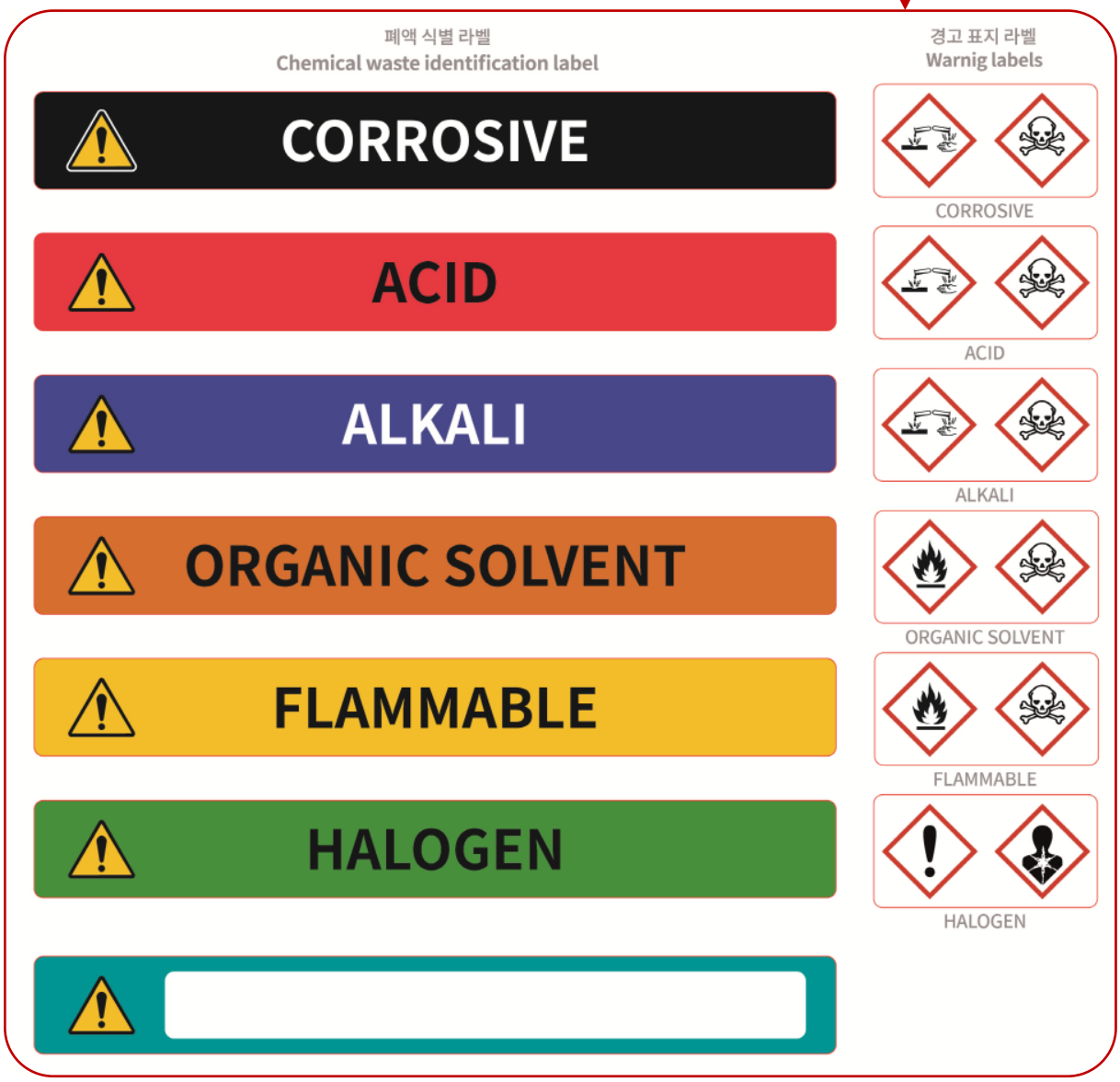
## 3.5. Other precautions

(1) Before collecting the explosive substances in the waste containers, please convert the reactive substances into safe substances.

(2) Please wear personal protective equipment when using and transporting the chemical waste containers.

(3) Store the chemical waste containers in a place with well-ventilated area and no direct sunlight.

(4) Store the waste container in a flat place so that it can be safely installed.

(5) Do not place it where people are busy, and dispose of waste containers as soon as possible in the laboratory.

**⚠WARNING**

• Do not expose the disposal system to the direct heat sources or direct sunlight.
• Store in a well-ventilated area.
• Keep the ambient temperature between 5°C and 40°C.
• The optimum temperature is 25°C.
• Dilute the acid or basic substances with water before collecting them.
• When collecting the original substances, you need to collect the secondary cleaning water in the waste container. However, when collecting highly toxic substances, you can collect the waste liquid with up to three times of cleaning water.
• Never place sediment or solids (gloves, bottles, tissues, magnets, chalk, etc.) in the waste container.

**⚠CAUTION**

• Place the waste container on a flat surface.
• Be careful not to place the waste container in a corridor where people often pass.
• Do not store more than three chemical waste containers(20L) in the laboratory.
• Do not exceed three months, the maximum storage period, for the waste container storage.
• When you handle the chemical substances by using or transporting, be sure to wear gloves, goggles, and other personal protective equipment in case of leakage.

©2025. All Rights Reserved. Jeio Tech Co., Ltd.

# 4. Specification

## 4.1 Safety funnels

| Model | SF-14C | SF-14P | SF-14Y | SF-14R | SF-14B |
|---|---|---|---|---|---|
| **Dimension (W x D x H, mm)** | 160 x 144 x 270 | | | | |
| **Opening diameter (ø, mm)** | 135 | | | | |
| **Cap size (ø, mm)** | 142 | | | | |
| **Volume (L)** | 1 | | | | |
| **Material** | Polypropylene (PP) | | | | |
| **Color** | Black (Electrostatic Conductive) | Persian green | Yellow | Red | Blue |
| **Cat. No.** | AAAM6001 | AAAM6002 | AAAM6003 | AAAM6004 | AAAM6005 |

## 4.2 Waste containers

| Models | WT-20C | WT-20P |
|---|---|---|
| **Dimension (W x D x H, mm)** | 228 x 300 x 402 | |
| **Opening diameter (ø, mm)** | 45 | |
| **Cap size (ø, mm)** | 78 | |
| **Volume (L)** | 22 | |
| **Material** | High-density Polyethylene (HDPE) | High-density Polyethylene (HDPE) |
| **Color** | Black (Electrostatic Conductive) | White |
| **Cat. No.** | AAAM6011 | AAAM6012 |

## 4.3 Accessories

| Description | Dimension | Cat. No. |
|---|---|---|
| **Carbon Filter** | 77 x 150 (ø x H, mm) | AAAM6051 |
| **Tray** | 370 x 470 x 100 (W x D x H, mm) | CBN0010 |
| **Earth plug set** | - | AAAF35453 |

# 5. Checking the product

(1) Please check the following components after opening the packaging of the product.

(2) If any damaged or missing parts or components are found, please contact the JEIOTECH sales department or seller.

(3) Safety funnel items

| Item | Quantity | Description |
|---|---|---|
| Safety funnel | 1 | Separate purchase available |
| Chemical waste identification label for safety funnel | 1 | Included with a safety funnel |
| 70% Level bar | 1 | Included with a safety funnel |
| Carbon filter & Filter timer | 1 | Optional, separate purchase available (replace every 6 months) |
| Earth cable with clamp | 1 | Included with safety funnel (Only available on SF-14C models) |
| Earth plug set | 1 | Optional, separate purchase available |
| User Manual | 1 | Included with the product |

(4) Chemical waste container items

| Item | Quantity | Description |
|---|---|---|
| Chemical waste container | 1 | Separate purchase available |
| Chemical waste identification label for chemical waste container | 1 | Included with a safety funnel |
| Chemical waste container label | 1 | Included with a chemical waste container |
| Tray | 1 | Optional, separate purchase available |
| User Manual | 1 | Included with the product |

©2025. All Rights Reserved. Jeio Tech Co., Ltd.

# 6. Warranty

## 6.1 Warranty period

(1) If the product fails during normal and proper use, the warranty period for manufacturing liability shall be two years from the date of delivery.

(2) If you request a repair, please check the items below and let us know so that you can get a quicker and more accurate repair.

- Date of purchase:
- Customer name / address / phone number / E-mail:
- Fault status
- Model
- Serial number

## 6.2 Technical Service Contact Points

Technical Services
Address : 153 (Youngsan-dong), Techno 2-ro Yuseong-gu, Daejeon-city, 34025, Republic of Korea
Tel : +82 (0)2-2627-3824
Website : www.jeiotech.com/eng/
E-mail : overseas@jeiotech.com

## 6.3 Certificate of Product Warranty

This is a guide to the product quality warranty and service terms provided by Jeio Tech Co., Ltd

All products are guaranteed to operate normally when installed under proper conditions, as specified in the user manual for each product, and used according to the intended purposes and permissible operating conditions.

The warranty period for product quality starts from the delivery date and is indicated in the catalog provided for each product or in the purchase order agreed upon with the customer. If a product malfunctions within the warranty period, Jeio Tech Co., Ltd. provides free repair services. However, among the transportation costs incurred during the repair process, any special costs resulting from the customer's site conditions (e.g., door or wall modifications, use of special handling equipment) are not covered under the warranty. Furthermore, Jeio Tech's liability under this warranty is limited to replacing the product or refunding the original purchase price at Jeio Tech's sole discretion if reasonable attempts to repair or replace the product fail.

Parts replaced during free repairs within the warranty period become the property of Jeio Tech Co., Ltd. Additionally, the warranty period is not extended due to free repairs. Jeio Tech is not responsible for any losses incurred by the customer due to delays in repairs caused by issues such as part shortages or shipping delays.

The following cases will be handled as paid repairs, even within the warranty period:

1. Non-defective cases

* Requests to resolve issues arising from improper operating environments or abnormal installations.

* Reinstallation requests or connections with external equipment other than the initial installation.

* Requests for usage instructions or simple adjustments that do not require disassembly.

* Inspections of non-defective products.

* Cleaning or debris removal service requests.

2. Failures caused by user errors:

* Malfunctions due to usage not in accordance with the user manual.

* Issues resulting from improper operating environments or abnormal installations.

* Malfunctions caused by incompatible chemical usage.

* Malfunctions caused by irregular or improper electrical or utility supplies.

* Damage due to product movement, storage, external impact, or dropping by the customer.

* Malfunctions resulting from unauthorized modifications or repairs without prior written approval from Jeio Tech Co., Ltd.

©2025. All Rights Reserved. Jeio Tech Co., Ltd.

* Malfunctions caused by using consumables or optional components not specified by Jeio Tech Co., Ltd.

* Repairs performed by personnel other than Jeio Tech or its authorized service agents.

3. Other cases

* When the product model or serial number has been damaged, removed, manipulated, or altered.

* Malfunctions caused by fire, flooding, or other natural disasters.

* Failure due to consumable parts reaching the end of their lifespan (e.g., lamps, filters, fuses, gaskets, packings, hoses, heat transfer fluids).

* Natural wear and tear (e.g., discoloration, screen printing degradation).

* Damage or malfunction of product accessories (e.g., power cables, data cables, glass components) or options.

* Malfunctions caused by errors in computers, computer components, or connected network systems.

• Before Jeio Tech performs free or paid repairs or accepts the return of a product, the customer must ensure the product has been decontaminated of chemical and biological hazards.

• Replacement parts used during repairs may differ from the original but will function correctly and be equivalent in performance.

• Jeio Tech does not bear any responsibility for damage to stored samples, data loss, or any other direct or indirect damages caused by product malfunctions, regardless of the warranty period.

• Jeio Tech does not guarantee issues related to computers, computer components, or network systems connected to the product, regardless of the warranty period.

• This warranty is the sole and exclusive warranty provided. Any implied warranties imposed by applicable laws are limited to the warranty period stated herein. This warranty cannot be modified by any verbal advice, statements, or documents on behalf of Jeio Tech unless explicitly labeled as a "Warranty Amendment" and signed by an authorized representative of Jeio Tech.

• This warranty applies only to products purchased for direct use by the customer. It does not cover products purchased for resale purposes without prior written approval from Jeio Tech.

• This warranty does not apply to consumables and is non-transferable.

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